



## Herasafe 2025 Biological Safety Cabinet

Smarter protection that never takes a day off

# Smarter protection that never takes a day off

The most important factor in selecting a biological safety cabinet is the assurance of excellent containment. Thermo Scientific™ Herasafe™ 2025 Biological Safety Cabinets (BSCs) feature advanced technologies for maximum safety and protection that never take a day off.

Robust construction and excellent ergonomic design are combined with autonomous safety features and a simple-to-understand interface that always gives a clear status indication. The Herasafe 2025 BSC delivers reliable performance with superior precision for a broad range of applications requiring assured safety for product, personnel and environment.

Thermo Scientific™ biological safety cabinets are designed to maximise sample protection and user safety, with an emphasis on containment, comfort and convenience, all day, every day. For more than five decades, we have been a pioneer and innovator of clean air technologies. Our global vision and experience have enabled us to continually provide improved application-based solutions ideally suited for your requirements and working environment, with compliance to all major worldwide safety standards.



## **Precise self-adjusting protection**

Autonomous safety features maintain proper air flow and monitor critical conditions in real time without user interaction. This supports valuable research by enabling our customers to focus where needed.



## **Simplified operation and ergonomics**

An intuitive, simple-to-use interface with clear icons provides the information you need, while the electronic sash with fully hinged opening enables easy access for supplies and equipment, and simplifies cleaning and disinfection.

# Smart technology offers optimal containment and energy efficiency

Autonomous safety features like the Thermo Scientific™ SmartFlow™ Plus and DAVE™ Plus technologies help ensure continuous accuracy that never takes a day off—providing optimal protection from contamination by compensating to match the way you work and your lab environment

Sash opening: 200 mm | Efficient design  
SmartFlow Plus | DAVE Plus Technologies



## **Electronic front sash**

With one-touch simplicity, the front window automatically positions itself to help ensure proper working conditions. The sloped front window allows a fatigue-free working position and has an aerosol-tight seal when closed.

## **200 mm sash opening**

Optimum balance of comfort and efficiency.

# Simplified operation and ergonomics



A readily understandable graphical display provides the information you need, while the electronic sash with fully hinged opening enables easy access for supplies and equipment, as well as effortless cleaning and disinfection.

Helping ensure that the Herasafe 2025BSC is managing safety conditions autonomously and is simple and intuitive to use allows our customers to fully focus on a broad range of demanding applications.

## Intuitive | Simple | Convenient

Simple icons showing status

Toggle between 4 fields of information

1. Hours of use  
2. Inflow velocity  
3. Downflow velocity  
4. Time

Simplified

Push-button operations

Window operation



### Thermo Scientific™ Smart Clean Plus Technology


Fully opening front-hinged window design enables easy cleaning and disinfection with full access to the interior cabinet space.

Full opening also allows access to the biological safety cabinet for thorough cleaning and visible reach to all cabinet surfaces including rear wall.



### 304 stainless steel interior

The 304 stainless steel interior provides resistance to commonly used chemicals and regular cleaning protocols. The finish is “brushed and sanded,” which both reduces glare and offers better scratch resistance.

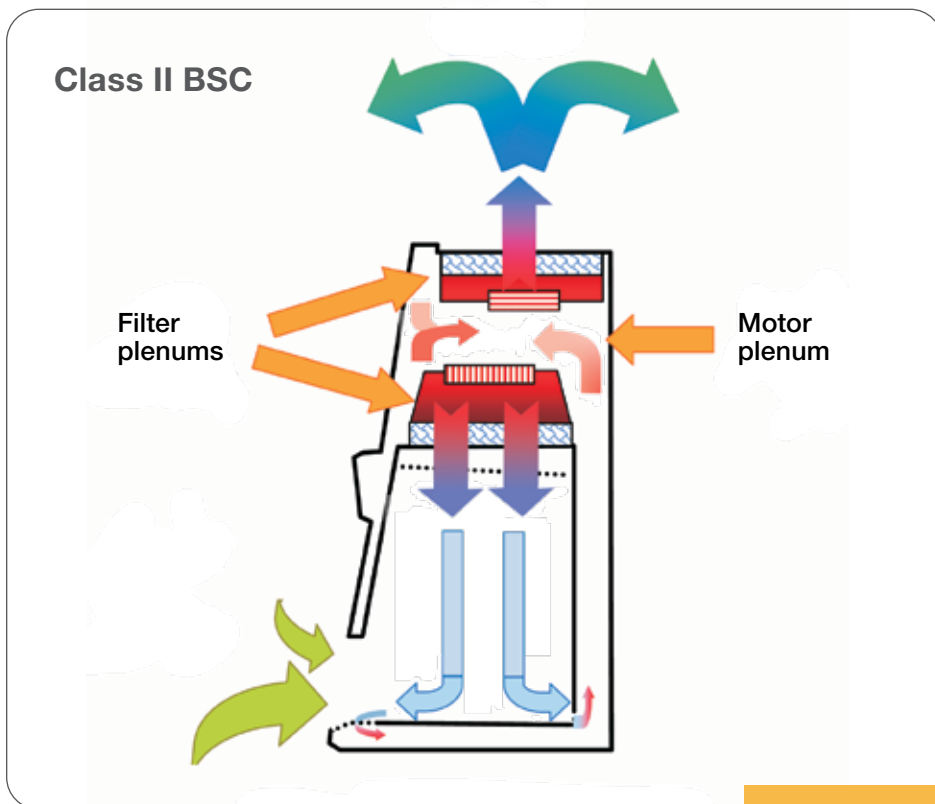

### SmartFlow Plus Technology

The SmartFlow Plus Technology automatically balances the performance of dual DC blower motors in order to maintain appropriate air flow rates whilst compensating for the loading condition of the protective HEPA filters unlike typical biological safety cabinets, which rarely perform at the same levels day after day due to HEPA loading, and potential disruptions to airflow such as paper covering vent openings.

#### How does this work?

Independent intelligent DC motors for downflow and inflow monitor and control fan speed in real time to maintain balanced airflow during filter loading or temporary obstruction.

Exclusive to our BSCs, independent supply and exhaust blowers automatically balance both the downflow and inflow velocities to maintain both operator and sample protection.

### DAVe Plus Technology

Our unique Digital Airflow Verification system validates protection by monitoring airflow velocities with independent pressure sensors that detect and alert to any changes across the exhaust and downflow plenums.

#### How does this work?

By sensing pressure in the exit of the blowers, DAVE Plus Technology monitors airflow directly for the inflow and downflow in the Herasafe 2025 BSC to help ensure both personnel and sample protection.



### Programmable crossbeam UV function

The Herasafe 2025 BSC is available with the option of an integrated powerful crossbeam UV disinfection system. The UV lights are installed on both side walls to illuminate the entire work surface without leaving shadows. This ensures thorough disinfection of the complete chamber for maximum sample protection.




### Energy efficient design

- Energy conserving dual DC motor design, generates less heat than older AC designs with double the usage lifetime and is digitally controllable for precision performance
- Stand-by mode automatically reduces blower speed to protect samples with reduced energy usage until being ready to resume work

### Optimising total cost of operation

SmartFlow Plus airflow compensation and balanced filters have led to safer BSCs that are easier to use and significantly less costly to operate.\*

#### The more you use it, the more you save!

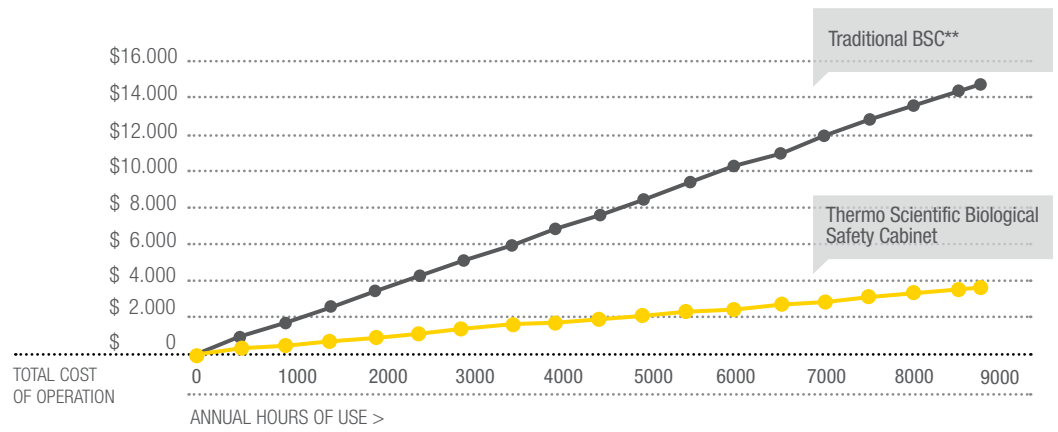

Costs, published in American Biotechnology Laboratory.

\*Compared to previous models.

\*\*Comparison based on 15 year life of 1.2 m wide biological safety cabinets.

"Traditional BSC" estimates based on published product specifications for readily available BSCs using AC motor technology.



# Innovative productivity enhancements simplify use, comfort and cleaning



## Customisation

The Herasafe 2025 BSC can be modified using our full customisation program. We have vast experience in designing custom adaptations to our safety cabinets. We can design from scratch or provide frequently requested solutions.

## Large interior work space/premium 304 non-corrosive stainless steel

Conveniently large interior workspace provides ample room for incorporation of instruments and process optimisation devices including robotic arms. We utilize a “brushed and sanded” finish to provide maximum scratch resistance and non-reflective properties.

## Interior LED lights

Visibility within the cabinet is enhanced with the use of bright, energy efficient lighting, reducing the potential for fatigue often experienced with use of standard fluorescent lighting.

## One-touch master button

Conveniently positioned master button opens and activates the Herasafe 2025 BSC simultaneously, correctly positioning the window and initiating blowers and lights.

## Smart Clean Plus hinged front window

Our exclusive fully opening front hinged window design enables easy cleaning and full access to the cabinet interior work space for thorough cleaning/disinfection and incorporation of ancillary work tools and equipment as may be desired.

## Padded armrests

Adjustable ergonomic armrests are provided with each unit to assure proper arm support to reduce strain during prolonged work periods while extending reach within the cabinet. Removable gel filled comfort pads are available as an option to enhance the user experience.



Smart Clean Plus hinged front window





Removable gel-filled comfort pads



# Outstanding comfort and ergonomic design

Featuring advanced designs for easy operation and maximum safety, the Herasafe 2025 Biological Safety Cabinet is recommended wherever sample and user protection is paramount.

- 10° sloped cabinet front—for exceptional comfort and ergonomics
- Bright workspace (side windows)—maximises light and visibility
- Unique internal armrest—positions hands in the centre of the work area for optimum protection
- Adjustable base stands—enables height to be set to suit environment and users
- Footrest accessory available for fitting to base stand
- Spacious work area—enables uncluttered work area for increased safety
- Motorised sash opening height of 200 mm for optimised usability and efficiency
- Standard battery backup to enable closure during power failure scenario

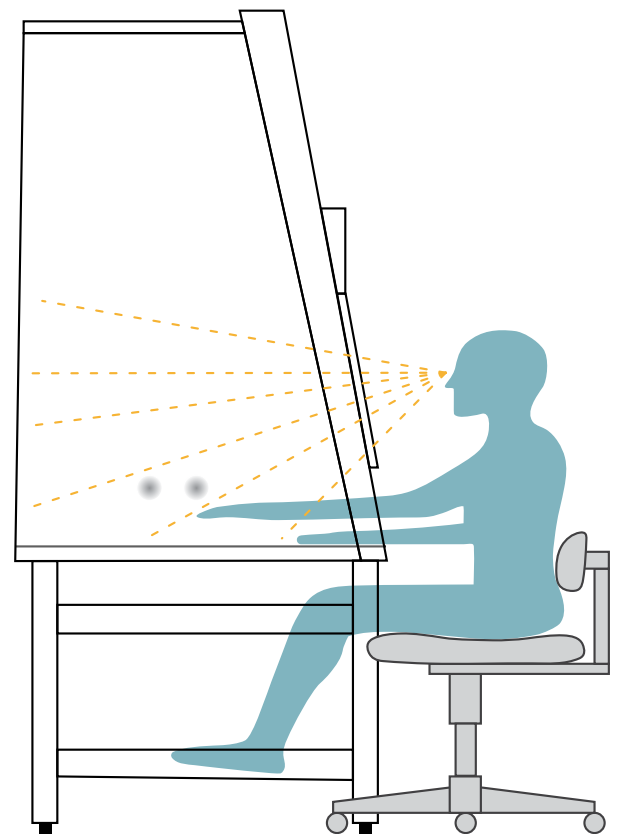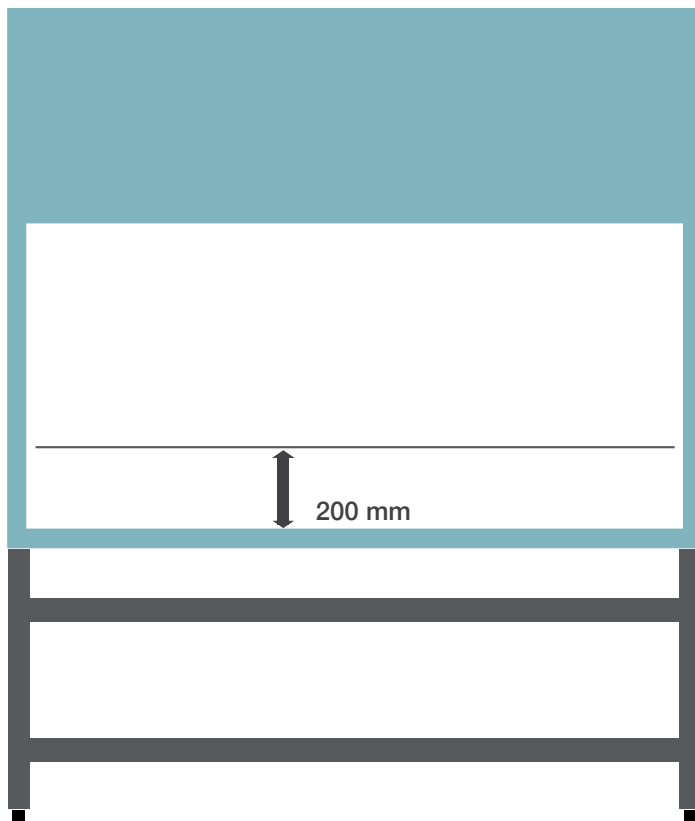



**Ordering information**
**Herasafe 2025 BSC**
**EN 12469**
**Dimensions**

Size/internal width in m (ft)			0.9 (3)	1.2 (4)	1.5 (5)	1.8 (6)
External dimension (cabinet only)	depth	mm (in.)	800 (31.5)	800 (31.5)	800 (31.5)	800 (31.5)
	width	mm (in.)	1000 (39.4)	1300 (51.2)	1600 (63.0)	1900 (74.8)
	height	mm (in.)	1536 (60.5)	1536 (60.5)	1536 (60.5)	1536 (60.5)
Usable work surface (net, without air intake grill)	depth	mm (in.)	465 (18.3)	465 (18.3)	465 (18.3)	465 (18.3)
	width	mm (in.)	900 (35.4)	1200 (47.2)	1500 (59.1)	1800 (70.9)
Total internal height of work area		mm (in.)	780 (30.7)	780 (30.7)	780 (30.7)	780 (30.7)
Shipping dimensions	depth	mm (in.)	925 (36.4)	925 (36.4)	925 (36.4)	925 (36.4)
	width	mm (in.)	1110 (43.7)	1410 (55.5)	1710 (67.3)	2010 (79.1)
	height	mm (in.)	1710 (67.3)	1710 (67.3)	1710 (67.3)	1710 (67.3)
Working height of front window		mm (in.)	200 (8)			
Maximum window opening (for cleaning, loading/unloading)		mm (in.)	773 (30.4)	773 (30.4)	773 (30.4)	773 (30.4)
Work surface height		mm (in.)	750-950 (29.5-37.4)			

**Weight**

Net weight	kg (lbs)	170 (375)	200 (441)	230 (507)	280 (617)
Shipping weight	kg (lbs)	190 (419)	225 (496)	260 (573)	315 (694)
Maximum load capacity of one-piece work surface	kg (lbs)	50 (110)	50 (110)	50 (110)	50 (110)
Maximum load capacity of single segmented work surfaces (per 30 cm segment)	kg (lbs)	25 (55)	25 (55)	25 (55)	25 (55)

**Electrical data**

Voltage	V	230			
Hertz	Hz	50/60			
Protection class		1/IP20			

**Filter specification**

Specification	H14 HEPA EN1822, 99.995% MPPS				
---------------	-------------------------------	--	--	--	--

**Certifications**

Standard	EN 12469				
Agency certified		TÜV	TÜV	TÜV	TÜV

**Ergonomics and utilities**

Lighting intensity	lux (fc)	1000 (93)	929 (86)	1500 (139)	1418 (132)
Interior outlets	pcs.	2 rear wall		2 rear wall	
Service valve ports	pcs.	2 on each side			
Slope of front window	degree	10			

**Herasafe 2025 noise, power, and heat specifications**

Noise level (measured according to EN 12469) 200 mm opening	dB (A)	57		57	
Power consumption (in normal operation, fans on, lights on) 200 mm opening	Wh	133	153	217	251
Power consumption (in stand-by mode)	Wh	55	55	80	80
Normal heat output in operation mode (non-vented) 200 mm opening	W	133	153	217	251
Stand-by heat output (lights off, non-vented)	W	55	55	80	80

**Exhaust volumes**

Exhaust volume of cabinet 200 mm opening	m <sup>3</sup> /h (cf/m)	288 (169)	386 (227)	483 (284)	580 (341)
Exhaust volume, thimble ducting 200 mm opening	m <sup>3</sup> /h (cf/m)	385 (226)	513 (302)	641 (377)	770 (453)

## Ordering information

Herasafe 2025 BSC	EN 12469			
Size/internal width in m (ft)	0.9 (3)	1.2 (4)	1.5 (5)	1.8 (6)
200 mm sash with cross beam UV-C in side walls	51033315	51033316	51033317	51033318
200 mm sash – no options included	51033311	51033312	51033313	51033314

### Receptacle options

Country-specific configurations are available to ensure correct line cord and internal receptacles are shipped with the cabinet. Standard version for EN12469 cabinets is Schuko, and for NSF cabinets is type B.

Switzerland	Option	51900300
UK	Option	51900303
Italy	Option	51900306
FR/B/CZ/SK/PL	Option	51900921
Australia	Option	51900449
Denmark	Option	51900481
Israel	Option	51901255

### Stand options

We recommend ordering a base stand to ensure a robust and stable installation of the BSC. All of our base stands have been tested with our units for stability and functionality. Choose between a fixed height, manually adjustable working height or electrically adjustable working height. Please consider the total height of the room when choosing the right working height. Footrest and castor options are available for added convenience.

Fixed height stand, charcoal grey Provides work height of 750 mm	Accessory	50155665	50155666	50155667	50155668
Adjustable height stand, charcoal grey Provides work height of 750-950 mm	Acc.	50155688	50155689	50155690	50155691
Electric adjustable stand, grey-white (RAL 9002)	Acc.	50117449	50109314	50117450	50109325
Footrest bar (for non-electric floorstand)	Acc.	50051982	50051983	50051984	50051985
Castor-set for stand, grey-white (4 pcs.) (RAL9002) Provides work height of 850 mm (for fixed stands only)	Acc.	50109977			

### Work surface options

You have the option to fully customise your new BSC to accommodate any work task. The standard EN cabinet includes our stainless steel segmented worktray for easy removal, autoclaving, and cleaning under the work surface. The standard NSF cabinet includes our stainless steel single piece indented worktray. To customise your work area, order the part number to remove the standard worktray and configure the cabinet as desired.

Without table tops (negative price)	Option	51900651	51900652	51900653	51900654
Flat work surface module (150 mm)	KAM15	Acc.	50073666		
Flat work surface module (300 mm)	KAM30	Acc.	50072685		
Flat work surface module (600 mm)	KAM60	Acc.	50073667		
Flat work surface module (900 mm)	KAM90	Acc.	50076398		
Flat work surface module (1200 mm)	KAM120	Acc.	50082970	50082970	
Flat work surface module (1500 mm)	KAM150	Acc.	No		50122210
Flat work surface module (1800 mm)	KAM180	Acc.	No		
Lowered work surface module (600 mm)	KAML60	Acc.	50073668		
Lowered work surface module (900 mm)	KAML90	Acc.	50073669		
Lowered work surface module (1200 mm)	KAML120	Acc.	No	50073670	
Lowered work surface module (1500 mm)	KAML150	Acc.	No		50073671
Lowered work surface module (1800 mm)	KAML180	Acc.	No		
Vibration dampened surface module with granite plate (300 mm)	KAMW30	Acc.	50073757		
Vibration dampened surface module with granite plate (600 mm)	KAMW60	Acc.	50076162		

## Ordering information

Herasafe 2025 BSC				EN 12469			
Size/internal width in m (ft)				0.9 (3)	1.2 (4)	1.5 (5)	1.8 (6)
Replacement armrests KAS 30 (Set of 2 armrests 316 SS)		KAS 30	Accessory	50073663			
Padded armrest 316 SS (Set of 2 including gel inserts)			Acc.	50156438			
Replacement gel pads for padded armrest (2 pcs.)			Acc.	50156437			

### Additional filter options

Additional exhaust HEPA filters (AEFs) can be fitted to the top of the BSC to provide 99.99999975% cumulative exhaust filtration. Please note that AEFs are not a requirement for recirculating BSCs, even in the UK, where BS5726 has been superseded by EN12469, since 2000. Additional exhaust carbon filters can be fitted on top of the BSC, but should only be used for traces of solvents being used as an adjunct to microbiological study.

Additional exhaust HEPA filter, retro-fit		KAEF	Acc.	50078857	50073772	50078856
Additional exhaust carbon filter, retro-fit		KAEFC	Acc.	50078783	50073773	50078784
Grid filter set, G1 acc. to EN 779 (5 pcs)			Accessory	50076846	50076921	

### UV options

Cross beam UV-C light option in side walls includes protective cover

Included in models: 51033315, 51033316, 51033317, 51033318

Standard

Standard (Cross Beam - in side walls)

Available for models: 51033311, 51033312, 51033313, 51033314

Option

51900642 (Cross Beam - in side walls)

### Media valves

Our standard BSCs feature a dual side-wall negative pressure design, which allows the user to install service valves through the side access ports. We also offer factory-installed piping options on our premium units to connect utilities to the top of the BSC to conserve right to left exterior footprint.

### Rear wall media connections

3 utility connectors prep (rear wall)

Select rear wall service taps additionally

Available for models: 51033311, 51033312, 51033313, 51033314

Option

51901251 (G3/8")

### Rear wall service taps/media valves (supplies are connected above BSC)

Combustible gas tap for rear wall		Acc.	50068235				
Non-combustible gas tap for rear wall		Acc.	50068234				
Water tap for rear wall		Acc.	50068233				
Vacuum tap for rear wall		Acc.	50068232				

### Side wall service taps/media valves/lead throughs

Combustible gas tap for side wall		Acc.	50046015				
Non-combustible gas tap for side wall		Acc.	50045959				
Water tap for side wall		Acc.	50044679				
Vacuum tap for side wall		Acc.	50044678				
Hose connection for side wall, hose dia. 10-13mm		Acc.	50059017				
Hose connection for side wall, hose dia. 3mm		Acc.	50076408				
Hose connection for side wall, hose dia. 6/3mm		Acc.	50066668				

### Solenoid valve (should be quoted with combustible gas G3/8 connection)

Solenoid valve for combustible gas tap		Acc.	50074905				
--	--	------	----------	--	--	--	--

### Exhaust duct connection kits

HEPA filters do a great job of trapping particles and biohazards, but they do not trap gases and toxic chemicals that may be used as an adjunct to microbiological studies. In these applications, it is important to select an exhaust duct transition to route exhaust air out of the laboratory to the external environment. We offer two different connections that can be installed in three different ways, depending upon facility requirements.

Thimble duct exhaust transition (Connection diameter 200 mm/8")		KTD	Acc.	50078796	50073770	50078795
Direct duct exhaust transition with anti blowback valve (Connection diameter 200 mm/8")		KDD	Acc.	50078791	50073771	50078790

## Ordering information

Herasafe 2025 BSC			EN 12469			
Size/internal width in m (ft)			0.9 (3)	1.2 (4)	1.5 (5)	1.8 (6)
User convenience options						
We offer unique customer-installed accessories to help organize your workspace and free up valuable space inside and around the Herasafe 2025 BSC.						
Hanging bar with set of 6 hooks	IVB 6 6	Acc.	50073943	No		
Hanging bar with set of 12 hooks	IVB 10 12	Acc.	No	50073944		
Utility basket to be used with hanging bar	KB	Acc.	50078900			
2 utility baskets incl. additional hanging bar	KBB	Acc.	50078901			

## Ordering Information

Herasafe 2025 BSC			Cat. No.	m (ft)
Additional topboxes for fumigation				
KSF frame fixed to top of each recirculating cabinet to accept the relevant topbox for fumigation. Brackets are fitted to side of frame to allow the use of speed clamps. Height = 50 mm.			KSF9	Acc.
			KSF12/15	Acc.
			KSF18	Acc.
KSNB topbox for use with hydrogen peroxide fumigation. Speed clamps allow for fixture to a KSF frame. The 35 mm hole will accept a Steris OR Bioquell outlet camlock fitting. Height = 110 mm with gasket.			KSNB9	Acc.
			KSNB-12/15	Acc.
			KSNB-18	Acc.
KS-EVBOX for use with formonia system and formaldehyde. Speed clamps allow it to be clamped to KSF frame, while the 2 threaded holes accept the ABEKP3 charcoal filters. This box can only be used if the formaldehyde is being neutralised with ammonia after fumigation. Height = 110 mm with gasket.			KS-EVBOX-9	Acc.
			KS-EVBOX-12/15	Acc.
			KS-EVBOX-18	Acc.
KSDB for use with formaldehyde fumigation only. Speed clamps allow for fixture to a KSP frame. This top box is supplied with a 3 m flexi trunking that allows formaldehyde to vent out into another exhausted application to purge fumes from cabinet. Side exit only. Height = 190 mm with gasket.			KS9DB160	Acc.
			KS12DB160	Acc.
			KS18DB160	Acc.
KSP(FD) for use with formaldehyde fumigation. Fixes to a KSF frame. Supplied with a 3 m flexi trunking that allows venting of fumes to an exhausted application. Top exit only. Height = 50 mm (frame).			KSP(FD)-9	Acc.
			KSP(FD)-12/15	Acc.
			KSP(FD)-12/15	Acc.
ABEKP3 charcoal filters screw into KS-EVBOX. Pack is 3 filters.			KSP(FD)-12/15	Acc.

Not all BSCs are manufactured to the same safety standards. Our design philosophy is grounded in compliance with the recognised worldwide safety standards, because your work cannot compromise on safety. Our biological safety cabinets are subjected to the following certifications:



### EN 12469

Thermo Scientific Class II biological safety cabinets are independently tested and certified to EN 12469 standard for Class II biological safety cabinets by TÜV Nord, the leading testing agency in Europe.

## Uw dealer en servicepartner in NL:



**Poly Temp Scientific BV**  
De Marne 211 Bolsward  
info@polytemp.nl  
0515-575105

Find out more at [thermofisher.com/bsc](https://thermofisher.com/bsc)

**ThermoFisher**  
SCIENTIFIC

This product is intended for General Laboratory Use. It is the customer's responsibility to ensure that the performance of the product is suitable for customer's specific use or application. © 2020 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified. **HERASAFE2025 0221**